

# LeMaitre® Embolectomy Catheters

## ORDERING INFORMATION



### Single Lumen Embolectomy Catheter

Highly rupture resistant balloon with strong bonds for excellent performance in even challenging situations. 6 years shelf life.

Balloon Diameter	Balloon Volume	Diameter	Color	Length	Model # [REF]
4.5 mm	0.05 ml	2F		40 cm	1601-24
4.5 mm	0.05 ml	2F		60 cm	1601-26
4.5 mm	0.05 ml	2F		80 cm	1601-28
8.0 mm	0.20 ml	3F		40 cm	1601-34
8.0 mm	0.20 ml	3F		80 cm	1601-38
10.5 mm	0.75 ml	4F		40 cm	1601-44
10.5 mm	0.75 ml	4F		80 cm	1601-48
13.0 mm	1.50 ml	5F		40 cm	1601-54
13.0 mm	1.50 ml	5F		80 cm	1601-58
13.5 mm	1.60 ml	6F		80 cm	1601-68
14.0 mm	1.75 ml	7F		80 cm	1601-78



### Over-the-Wire Embolectomy Catheter

Highly rupture resistant balloon. Compatible with standard guidewires. Facilitates dye enhancement and irrigation option (e.g. heparin saline solution, thrombolysis, etc.). Stopcock and syringe included. 6 years shelflife.

Balloon Diameter	Balloon Volume	Diameter	Color	Length	Guidewire	Model # [REF]
6 mm	0.20 ml	3F		40 cm	0.018"	1651-34
6 mm	0.20 ml	3F		80 cm	0.018"	1651-38
10 mm	0.75 ml	4F		40 cm	0.025"	1651-44
10 mm	0.75 ml	4F		80 cm	0.025"	1651-48
12 mm	1.50 ml	5.5F		40 cm	0.035"	1651-84
12 mm	1.50 ml	5.5F		80 cm	0.035"	1651-88
13 mm	1.60 ml	6F		40 cm	0.035"	1651-64
13 mm	1.60 ml	6F		80 cm	0.035"	1651-68
14 mm	1.75 ml	7F		80 cm	0.038"	1651-78



### NovaSil™ Silicone Single Lumen Embolectomy Catheter\*

Due to the Silicone balloon material the Novasil catheter provides a latex-free option for latex sensitivity. 5 years shelflife.

Balloon Diameter	Balloon Volume	Diameter	Color	Length	Model # [REF]
4 mm	0.05 ml	2F		40 cm	1801-24
4 mm	0.05 ml	2F		60 cm	1801-26
6 mm	0.20 ml	3F		40 cm	1801-34
6 mm	0.20 ml	3F		80 cm	1801-38
9 mm	0.60 ml	4F		40 cm	1801-44
9 mm	0.60 ml	4F		80 cm	1801-48
11 mm	1.00 ml	5F		40 cm	1801-54
11 mm	1.00 ml	5F		80 cm	1801-58
13 mm	1.60 ml	6F		80 cm	1801-68
14 mm	1.75 ml	7F		80 cm	1801-78

\*Not made with natural rubber latex



Your Peripheral Vision™

LeMaitre is a registered trademark and NovaSil and Your Peripheral Vision are trademarks of LeMaitre Vascular, Inc.  
©2013 LeMaitre Vascular, Inc. All rights reserved. M0005 Rev. J 06/13



### LeMaitre Vascular GmbH

Otto-Volger-Str. 5a/b  
65843 Sulzbach/Ts.  
Germany

T +49 (0)6196 659230  
F +49 (0)6196 527072

### LeMaitre Vascular SAS

79 avenue de Villiers  
75017 Paris  
France

T +33 (0)3 44 26 00 41  
F +33 (0)3 44 26 07 31

### LeMaitre Vascular S.r.l.

Via Leone Tolstoj 86  
20098 San Giuliano Milanese (MI)  
Italy

T +39 02 988 48 51  
F +39 02 984 97 059

### LeMaitre Vascular Spain S.L.

Calle Basauri 17  
28023 Madrid  
Spain

T +34 902 500 037  
F +34 91 184 98 71

### LeMaitre Vascular Switzerland GmbH

Grabenstrasse 25  
CH-6340 Baar  
Switzerland

T 0800 561 761  
F +41 415 608 236

### LeMaitre Vascular GK

1F Kubodera Twin Tower Bldg.  
2-9-4 Kudan-minami, Chiyoda-ku  
Tokyo 102-0074  
Japan

T +81 (0)3 5215 5681  
F +81 (0)3 5215 5682

### LeMaitre Vascular ULC

5915 Airport Road, Suite 608  
Mississauga, Ontario  
Canada L4V 1T1

T +1 905 673 2266 0  
F +1 905 673 2223

### LeMaitre Vascular, Inc.

63 Second Avenue  
Burlington, MA 01803  
USA

T +1 781 221 2266  
F +1 781 221 2223



www.lemaitre.com

# LeMaitre<sup>®</sup> Embolectomy Catheters

The World's Strongest Embolectomy Catheter



 **LeMaitre<sup>®</sup>**  
VASCULAR

Your Peripheral Vision<sup>™</sup>

# LeMaitre<sup>®</sup> Embolectomy Catheters

**STRONGER BOND. STRONGER BALLOON.  
BEST VALUE.\***

#### Superior strength.

Significantly stronger Single Lumen Embolectomy Catheter than leading competitor\*. It features stronger bonds, which allow to apply more force with less risk of balloon detachment from the catheter and a stronger balloon to yield greater burst strength and reducing the risk of balloon rupture in even challenging situations.

#### Extended Shelf Life

LeMaitre Embolectomy Catheters have now up to 6 years shelf life.

#### Trusted brand

LeMaitre Embolectomy Catheters have been sold internationally for over 20 years and have been used in over 1 million procedures.

#### OVER-THE-WIRE EMBOLECTOMY CATHETER

##### Guidance.

Compatible with standard guidewire to direct catheter to proper position

##### Dye Enhancement.

Inject contrast media to facilitate emboli/thrombi imaging

##### Irrigation.

Allows for lubrication during catheter insertion and withdrawal



#### LONGER SHELF LIFE

Latex: 6 yrs  
Silicone: 5 yrs

## STRONGER BOND TEST RESULTS\*

N = Force in Newton

Mean Bond Strength in N +/- Standard Deviation

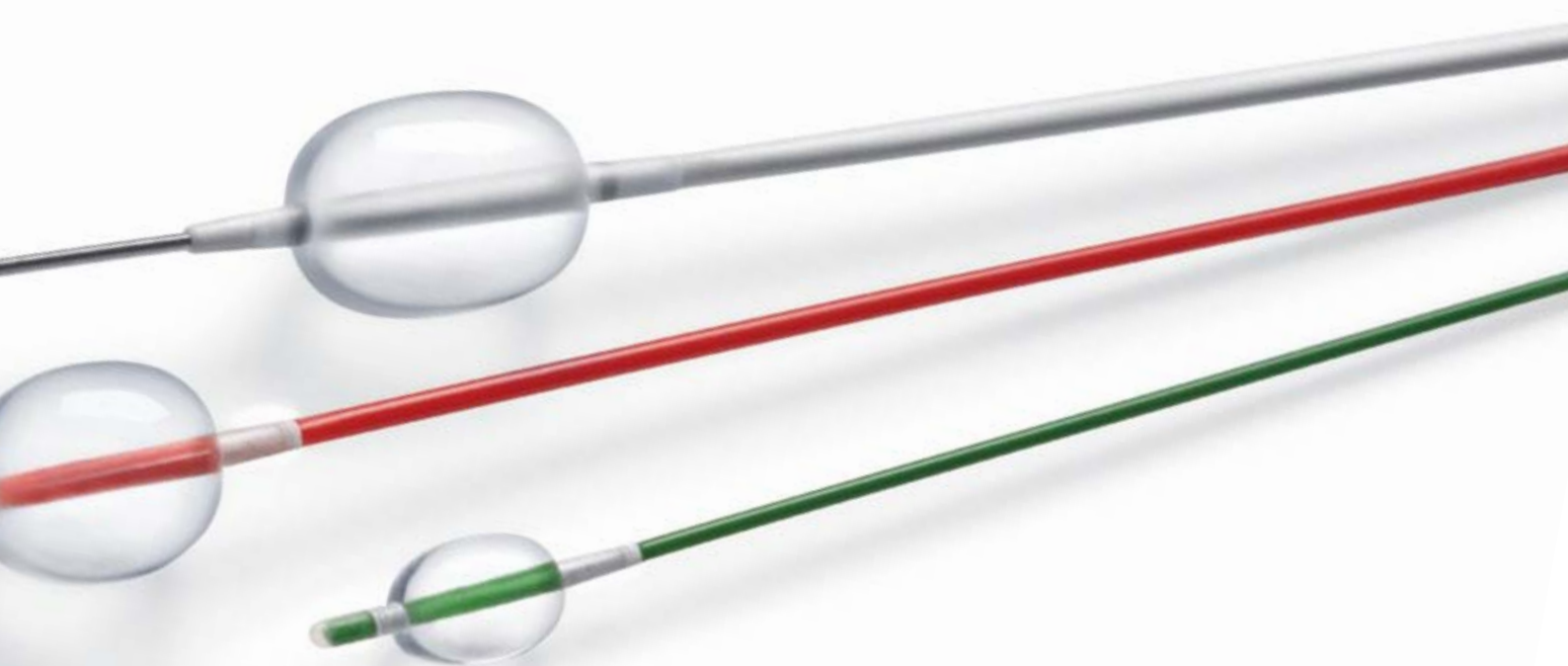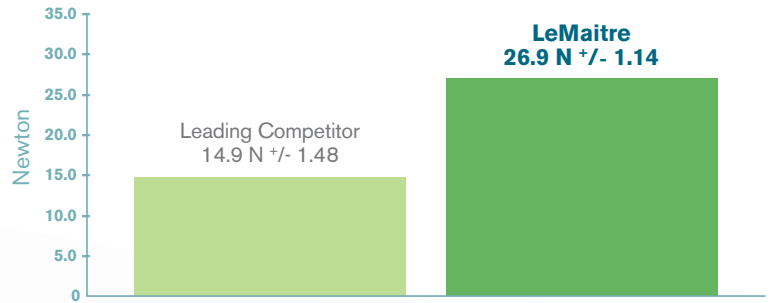

### SINGLE LUMEN EMBOLECTOMY CATHETER

60% stronger, rupture-resistant balloons\*  
80% stronger catheter bonds\*

Stronger bonds and balloons allow  
you to apply more force with less risk.

### NOVASIL SILICONE SINGLE LUMEN EMBOLECTOMY CATHETER\*\*

The NovaSil Silicone balloon  
catheter provides a viable option  
for patients with latex sensitivity.

\*Comparative bench top testing performed on the LeMaitre 4F Single Lumen Latex Embolectomy Catheter and the leading competitor's 4F catheter demonstrated up to 80% stronger balloon bonds and up to 60% stronger balloons.

\*\*Not made with natural rubber latex